


Niño TWBNL (Lateral) Series

TWBNL (Lateral)




Universal Track Mount

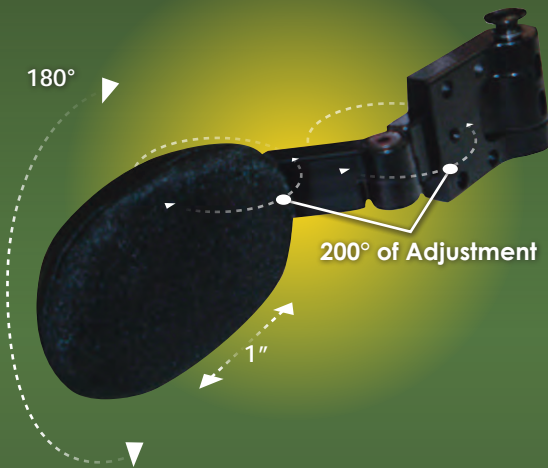

Can Mount
sizes: 1" - 7/8" - 3/4"



Stealth Track

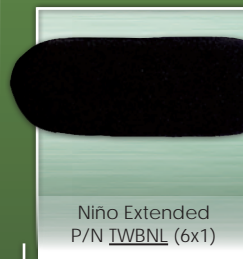

Niño Extended
P/N TWBNL (6x1)




Niño Standard
P/N TWBNL (3x3)



Niño Small
P/N TWBNL (3x2)

The Dynamic Full Surface Contact (FSC) Technology

Full surface contact is a desirable goal for Seating and Positioning. Often we are dealing with asymmetrical clients and symmetrical products. This inherent mismatch compromises the goal every time. Now, Stealth has developed seating and positioning components that have the ability to adjust to the body's real needs so contact between the component and the client is made full surface.

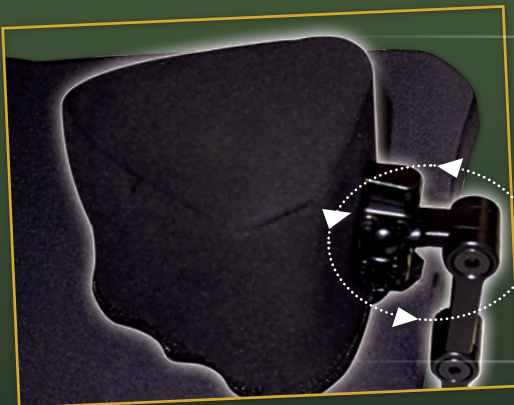

ABD-FSC ABductor

Once again, typical abductors are symmetrical devices and many clients are asymmetrical and need specifically angled pads/pommels. Stealth has solved this by developing the ABD-FSC ABductor. The ABductor has infinite angle capability and flip down hardware.



The ASTWBL Laterals for "All Seasons"


New to Stealth is the dynamic width adjustable Lateral. Needing no tools, very quick and totally intuitive. The ASTWBL Laterals are adjustable as and when needed, infinitely over a 4" range per side (8" overall per two laterals). The ASTWBL's can be used with the regular TWBL or the New TWBL-FSC's.

Niño i2i Series


The i2i Head and Neck Positioning for the most Challenging Clients!

The i2i was designed to be used as a therapeutic tool to encourage active head control within a limited range. It offers integrated support to the entire cervical spine (posteriorly, laterally and anteriorly, if the chin prompt is needed), the sub-occipital area (to provide a resting spot for the head), and the anterior upper thoracic area, which often requires cueing and support to promote good alignment. It has been successful with clients at either end of the spectrum:

- Those dominated by increased flexion, or too weak to hold their head up for any sustained amount of time.
- Those presenting, with a great deal of extensor tone but poor midline control, which often results in their head getting caught under or to the side of their existing system.

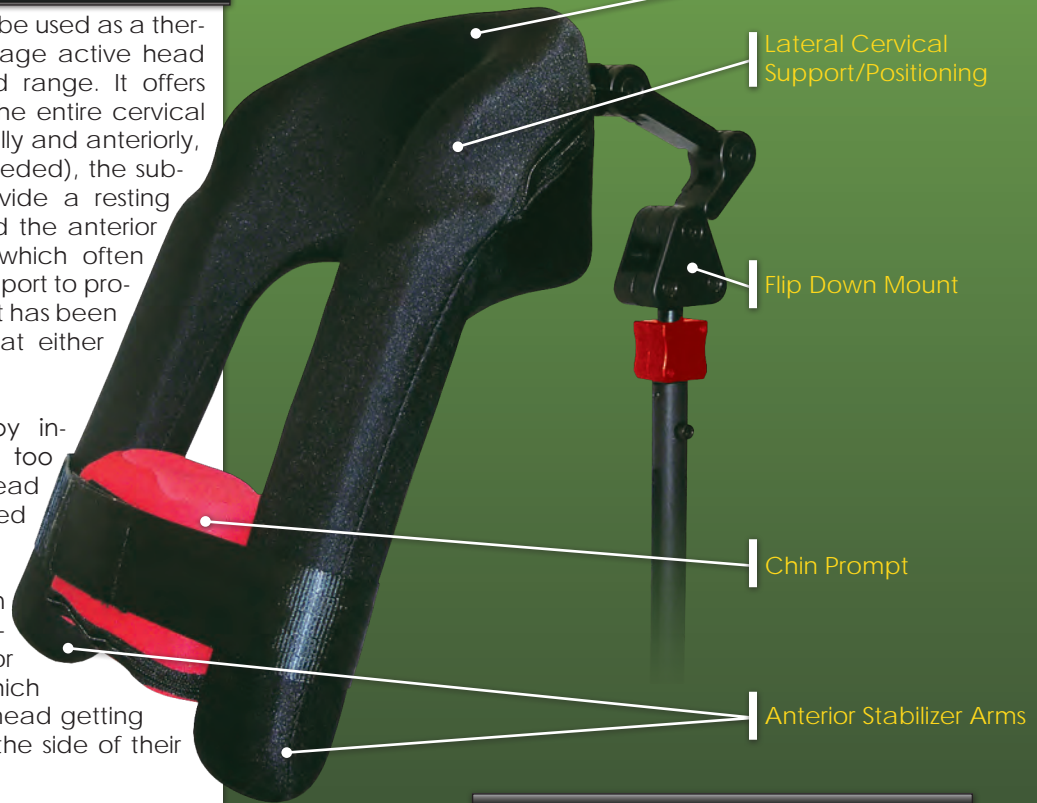

Removable Meal Time covers in vinyl are now provided with all i2i chin prompts.

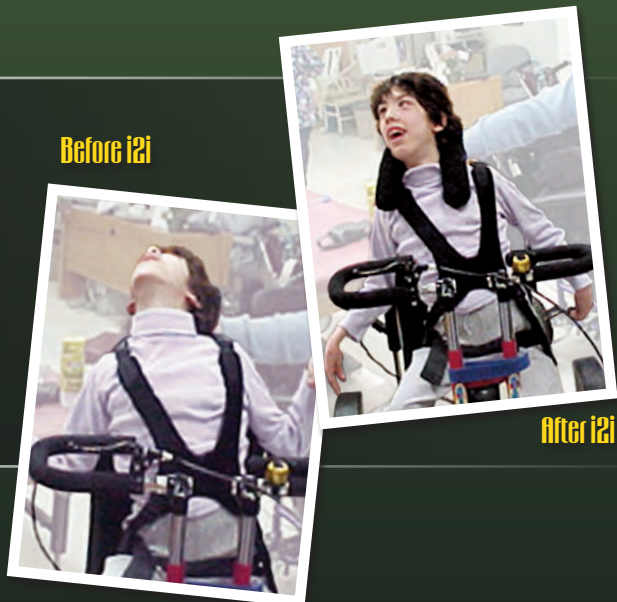

Niño i2i Series



Part Numbers	Maximum Neck Circumference	Overall Width	Arm Length	Maximum Neck Width
i2i 480 Small Standard	11"	8"	6.5"	4"
i2i 481 Small Short Arm	11"	8"	4.5"	4"
i2i 580 Medium Standard	14"	9"	7.5"	4.5"
i2i 581 Medium Short Arm	14"	9"	5.5"	4.5"
i2i 490 Small Low Profile	11"	8"	6.5"	4"
i2i 491 Small Low Profile Short Arm	11"	8"	4.5"	4"
i2i 590 Medium Low Profile	14"	9"	7.5"	4.5"
i2i 591 Medium Low Profile Short Arm	14"	9"	5.5"	4.5"

Fabric: (i2i headrest) Dartex, reverse Dartex, Neoprene, and natural sheepskin; (Chin Prompt) Neoprene or vinyl.

Comfort Plus Series

Niño
Comfort Plus
with HMS450
Pad Size: (6" L x 2.75" H)



Assembly Disengaged

CPN360 with SUNS9-L/E8675/ESYEL
and 910 Pad

Lateral Swing-Away Assembly

Works with Niño Ultra and CPN Comfort Series. A wide range of accessories and options are available.

SUNS9-B = Both Sides
SUNS9-R = Right Side Only
SUNS9-L = Left Side Only
*(Standard Facial Pads for Niño)

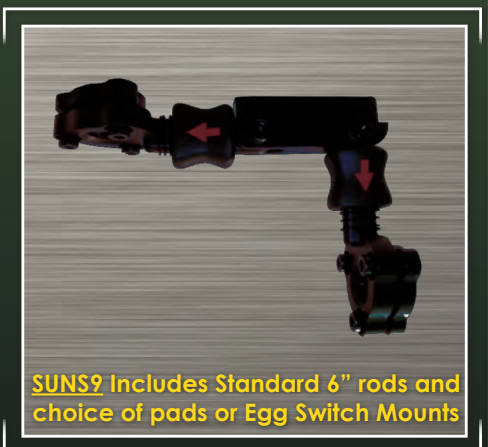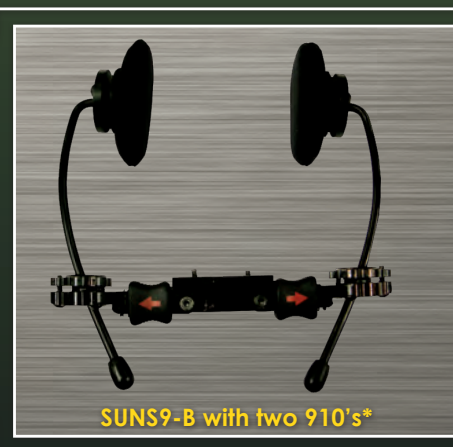

Niño Ultra QCR Series

The Niño Ultra QCR Series is an all encompassing Head Positioning System with adjustment capabilities and a plethora of options and accessories to meet virtually any client's need.



Sub-Occipital and Occipital incorporates independent fore and aft adjustments. Simple, obvious and reliable. Adjustment range is 3" with the Multi-Axis Mount and 5+" with the TWB LINK Mount

In addition, both Occipital and Sub-Occipital pads incorporate rotational, posterior, anterior and oblique adjustments.



SUN120 with SUNS9-B/910 PADS



The ULTRA QCR is like having multiple hands working to position your client's head.

Niño Ultra QCR Series

Quick Set Up Simple and Intuitive
Clean Design See the person, not the hardware
Reliable Holds your settings, **No Excuses!**



SUN120 with 910 Pads
 SUNS9-B and HMS450




Niño Ultra with 920 Pad

Niño Ultra with 910 Pads

Niño Ultra Pads

530 Pad (Occipital)	610 Pad (Sub-Occipital)	910 Pad (Facial)	900 Pad (Facial)	920 Pad (Facial) (Bendable)	930 Pad (Facial) (Bendable)	ESM Egg Switch Mount
-------------------------------	-----------------------------------	----------------------------	----------------------------	--	--	--------------------------------

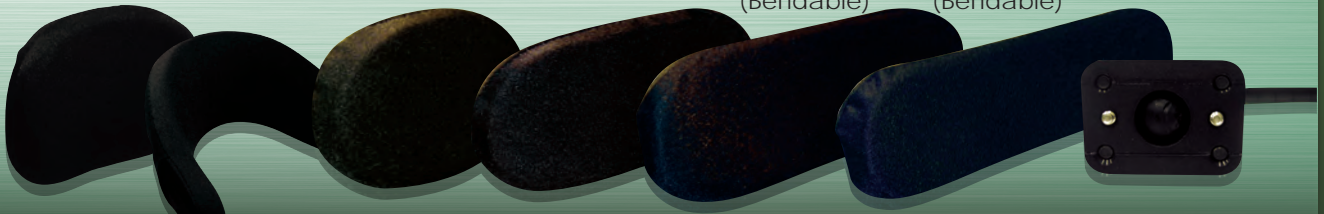

3" x 2"	6.375" L x 2.25" H	3" L x 2" H	4.5" L x 2" H	6.375" L x 2.25" H	6.375" L x 1.5" H	1.75" L x 1.25" H
---------	--------------------	-------------	---------------	--------------------	-------------------	-------------------

Lateral Swing-Away Assembly

A wide range of accessories and options are available.

SUNS9-B Both Sides
SUNS9-R Right Side Only
SUNS9-L Left Side Only

*(Standard Facial Pads for Niño)



SUNS9-B with two 910's*



SUNS9 Includes Standard 6" rods and choice of pads or Egg Switch Mounts

Niño SUNS9 Swing-Away



As Shown:
Swing-Away ([SUNS9-B] Right & Left) with a Long Lateral Pad (930) on the right and a Green Egg Switch (ESGRN) on the left.

Swing-Away arms offer a sleek, unobtrusive design while featuring an easy, one-hand release mechanism for swinging pads and switches back a full 100 degrees. This ensures convenience and is important for transfers, therapy and fatigue issues. Swing-Away pads or switches can be quickly and easily locked back into their original setting.

Facial Lateral Rod	
	
SU Rods	
SU260	Standard 9"
SU262	Standard for Niño Short 6"
SU264	13"
SU265	18"








F8675 (Figure 8) mounts are uniquely designed to provide 240 degree rotation capability for a second pad, switch or embedded speaker on a traditional Swing-Away.



SUNS9-B
Includes Standard 6" rods and choice of pads or Egg Switch Mounts

To Order Call:
1-800-965-9229
and for more info log on to:
www.stealthproducts.com

Facial Lateral Options

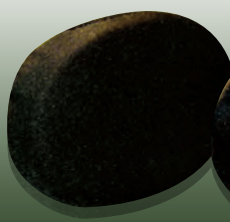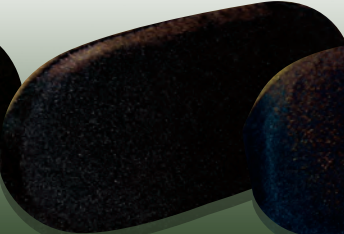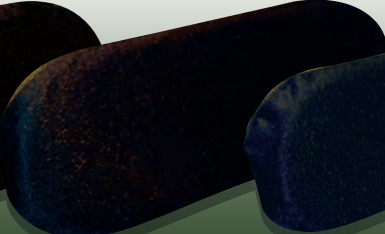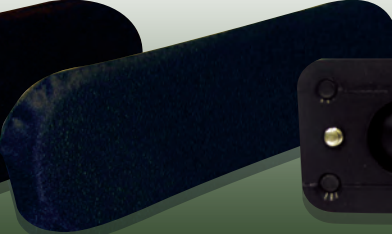910 Pad	900 Pad	920 Pad (Bendable)	930 Pad (Bendable)	ESM Egg Switch Mount
				
3" L x 2" H	4.5" L x 2" H	6.375" L x 2.25" H	6.375" L x 1.5" H	1.75" L x 1.25" H

Swing-Away Options

			
SM600 Mini proportional joystick mount	SM630 ASL MEC joystick mount	FO660 20" Flex fiber optic mount - Cable Included -	MC210 Tash min-cup switch
			 ML210 Micro Lite switch

			
FO650 Lateral rod fiber optic mount (Cable included)	CN600 12" Flex micro-lite switch mount	SP625 Sip & Puff mount	EGG SWITCHES
	does not include ML210		ESYEL ESRED ESBLK ESGRN ESBLU
			Also available: ES-EXT (Egg switch extension cable)