

All Positioning Series Headrests

TECHNICAL DATA

INCLUDES:

All Positioning Headrest Pad
17mm Mounting Ring

AP250

Small AP Headrest Pad

HMO-HB

Horizontal Bar

HMO475-TR

HMO475 Triangle


HMO475

Headrest Mounting Hardware

HMO475-VB

Vertical Bar

DIMENSIONS & MOUNTING

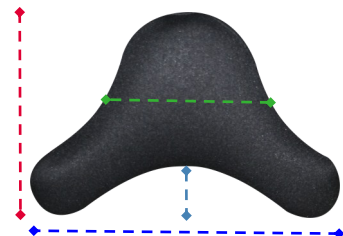

	Size	Width	Height		
AP150	X-Small	7"	6"	1 1/2"	5"
AP250	Small	8"	6 1/4"	3/4"	5"
AP350	Medium	9"	6 1/2"	1"	5"
AP450	Large	10"	7 1/2"	1 1/4"	6"
AP550	X-Large	13"	8"	1 1/2"	8 1/2"

17MM

BALL MOUNT

SDRING

17mm Mounting

17M BALL

MOUNTING:

The All Positioning headrest attaches to a 17mm ball mount via a SD-Ring. When you have the mounting ring in place, screw it into the adjacent holes on the headrest.

PRODUCT DESCRIPTION

When you need comfort and adjustability, depend on our All Positioning headrest. The AP Series offers a unique flexibility. The innovative memory metal allows an "adjust on demand" quality unlike most headrests. With five sizes available, all spectrums of shapes can be met.

TECHNICAL DATA

1. 5/32 tool for socket-head and flat-head screws on mounting ring and mounting hardware.
2. Includes: ball mounting ring and headrest pad.
3. Headrest attaches to: HMO400, HMO475, TWB480-1, TWB480-FDM, and TWB485-1.

SELLING POINTS

1. The AP's memory metal allows unique flexibility and adjustability.
2. Compatibility with all of our standard mounting hardware.
3. Ball mount allows for angle and rotational adjustment.

BUILD YOUR ALL POSITIONING HEADREST

STEP 1: PAD SIZE

The AP style headrest comes with several types of foam and covers. According to the individual's needs, determine the necessary pad size, whether foam or gel would be beneficial, and cover material. (See Headrest Pads on back)

AP150 — AP250 — AP350 — AP450 — AP550

FOAM/GEL

Open Cell Foam
Closed Cell Foam
Action Gel

COVERS

Dartex
Reverse Dartex
Vinyl
Sheep Skin
Neoprene

STEP 2: SWING-AWAYS (SUS9 SERIES)

Facial lateral Swing-Aways for pads, switches, speaker sites and chin control devices can be added to the AP headrest. (See Swing-Away Options on back)

SUS9-B-HW — SUS9DL-DB — SUNS9-B-HW

STEP 3: HARDWARE

The final step is deciding what type of hardware to use for mounting the AP headrest to your chair. (See Hardware Options on back)

HARDWARE

HMO400
HMO475
TWB480
TWB480-FDM
TWB485

HARDWARE OPTIONS



HMO400-17

The HMO400 is our straight mounting hardware can be adjusted fore and aft.



HMO475-17

The HMO475 allows fore and aft adjustment as well as rotational adjustment



TWB480-17MM

The TWB480 has a rotating triangle and links that can be adjusted to any position.




TWB485-17MM

The TWB485 combines the links of the TWB480 and curved bar of the HMO475.



TWB480-17FDM

An easy transfer flip down mount and links make Stealth's TWB480-FDM.

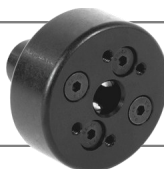

TD100

The Tone Deflector is design to absorb quick impact and avert slanted angle loading.

HEADREST PADS



AP150

Extra Small
All Positioning Headrest



AP250

Small
All Positioning Headrest



AP350

Medium
All Positioning Headrest



AP450

Large
All Positioning Headrest



AP550

Extra Large
All Positioning Headrest

SWING-AWAY OPTIONS

More information available at www.stealthproducts.com



SUS9-B-HW

Stealth Product's Swing-Away hardware offers a sleek, unobtrusive design while featuring an easy, one-hand release mechanism for swinging pads and switches back a full 100°.



SUS9DL-DB

The dual swing away design allows for more attachment sites for facial pads or switches for those clients who need more real estate or extra support.



SUNS9-B-HW

Convenient sizing, versatility, comfort and intuitive function with a clean, elegant design makes Stealth's SUNS9 the best choice when in need for pediatric facial pads and switches.



900
Facial Pad
(4.5 " x 2 ")



F8675

Stealth's figure 8 bracket adds real estate for an extra attachment.



910
Facial Spot Pad
(3 " x 2 ")



JOYSTICK MOUNTS

SM600: Mini Proportional mount
SM610: 1812 mount
SM620: Tash Mini mount
SM630: ASL MEC mount



920
Long Facial pad
(6 " x 2 ")



SP625

Sip & Puff mount



930
Long Niño Facial Pad
(6 " x 1.5 ")



CN600

Swing-away flex Micro-Lite switch.