


TWBNL—The World's Best Niño Laterals

TECHNICAL DATA

*Pad not included w/ part # TWBNL

TWBNL-L

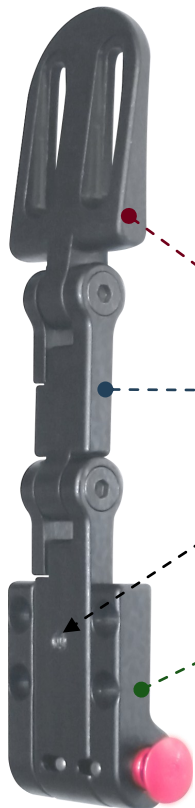

STANDARD CONFIGURATION

The TWBNL Includes:

- Pad Mount
- Multi Adjustable Link
- Quick Release Swing Away Mechanism.

PAD MOUNT

Niño Lateral Pad Mount—TWBNL-31

LINKS

Niño Multi-Adjustable Links—TWBNL-29

SET SCREW

Adjust set screw to prevent movement in the swing away mechanism

SWING-AWAY

Quick Release Swing Away Mechanism

QUICK RELEASE BUTTON

Spring Loaded Quick Release Button

TWBNL-R



TWBNL

Lateral Thoracic Articulating Adjustable Swing Away Hardware (E1028)

The Worlds Best Niño Laterals are an extremely versatile swing away pediatric thoracic lateral bracket that is designed to be fully adjustable to accommodate the needs of almost any client. The standard configuration comes with one quick release swing away mechanism, one extra link and one pad mount. A variety of different mounting hardware and multi-adjustable links are available to get any configuration you need.

PRODUCT DESCRIPTION

The ideal laterals need to imitate the movement of hands, arms and wrists so they can get to the 'right' spot and stay there! Just as the standard TWBL, the TWBNL is designed to do this as well! TWBNL Laterals 'snug' more easily and still retain strength and function. The World's Best Niño Laterals is Stealth's most versatile and adjustable pediatric solution to lateral needs.

TECHNICAL DATA


- $\frac{3}{32}$ " tool for $\frac{8}{32} \times \frac{1}{4}$ button head screws (pad mounting.)
- $\frac{1}{8}$ " tool for 10-32 $\times \frac{3}{16}$ " set screws and 10-32 $\times \frac{3}{4}$ " flat head screws (link joints).
- $\frac{5}{32}$ " tool for 10-32 $\times \frac{3}{4}$ socket head screws (hardware mounting.)

SELLING POINTS


- Push button quick release swing-away mechanism provides easy transfer from user to chair.
- Smooth finished edging on every piece.
- Fully Adjustable to accommodate the needs of almost any client.

DIMENSIONS

Side View

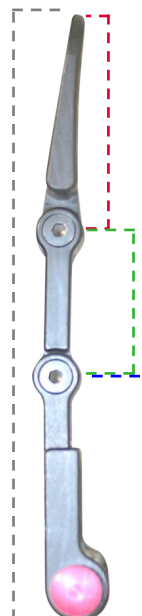

Front View

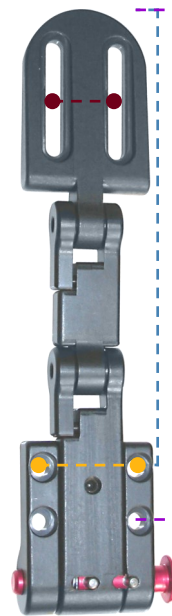


Standard Configuration Dimensions

Use These standard dimensions to determine which configuration you will need. See back for alternate mounting hardware and additional options.

Pad Mount	2½"
Link	1¾"
Swing Away Mechanism	2¾"
Pad Mount Slot Spacing	¾"
Swing Away Mechanism Hole Spacing	1⅞"
Mounting Hardware to Pad Mount	5⅝"
Mounting Hardware to Pad Mount	6"
Total Hardware Length	7"

2GTR Track Nut	LTNHW503
Niño Lateral Planar Mount	LTNHW505
All Season Bracket	LTNHW504

STANDARD HARDWARE ASSEMBLY

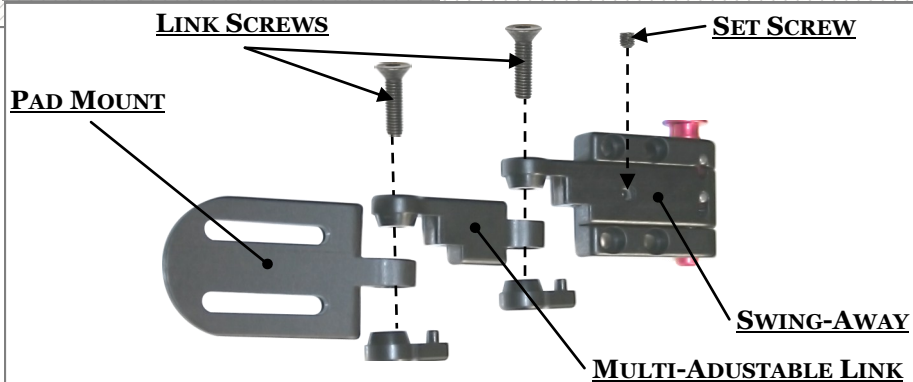

Example Configuration

TWBNL

This example configuration to the right is shown with an extra link and planar mount connected to a 7" track by anti-rotating nut.



Swing-Away Action

Stealth's adductors can be swung away 180° and helps in transferring the user from their chair.




More information available at www.stealthproducts.com

MOUNTING HARDWARE

Niño Mounting Hardware

Stealth offers mounting hardware of various sizes and style. Our diverse mounting hardware provides multiple options to mount on virtually any mobility base.

CANE MOUNT CLAMPS

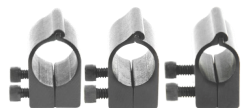

- 1" Niño Cane Clamp—LTNHW500
- 3/4" Niño Cane Clamp—LNTNW501
- 7/8" Niño Cane Clamp—LNTNW502

PLANAR MOUNTS*



- Niño All Season Bracket LTNHW503
- Niño Planar Mount LTNHW505

TRACKS



- 2" Track—SBHD-234
- 3" Track—SBHD-235
- 6" Track—SBHD-233
- 7" Track—SBHD-231
- 8" Track—SBHD-236
- 10" Track—SBHD-140
- 12" Track—SBHD-232
- 14" Track—SBHD-237
- 16" Track—SBHD-238
- 18" Track—SBHD-239



- Niño 2nd Generation Track Nut LTNHW503

INDIVIDUAL PARTS

Niño Links

Links can be added to increase hardware adjustability and length.

Niño Link Assembly

LNHW120




Assemblies

If you're ever in need of a simple replacement of a part, we have you covered!



Niño Left Swing Away Assembly

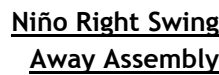

Niño Right Swing Away Assembly



Niño Lateral Pad Mount

LATERAL PADS



Niño Lateral Pads

With Stealth Product's thoracic positioning/support pads the sky is the limit. We not only offer a wide selection of standard size pads made from an aluminum base with 1" closed cell foam. We can also make any other size one might dream up with foams and covers of choice.†

† When ordering different upholstery, foam or non-standard size the pad is custom.

PADS



TWBNL2x3

2"x3" Niño Lateral Pad



TWBNL3x5

3"x5" Niño Lateral Pad



TWBNL2x3

6"x1" Niño Lateral Pad



TWBNL3x3

3"x3" Niño Lateral Pad