

TECHNICAL DATA

*Pad not included w/ part # TWBNL

STANDARD CONFIGURATION


The TWBTL Includes:

Pad Mount

Thin Profile Quick Release Swing Away Mechanism

TWBTL-L

TWBTL-R

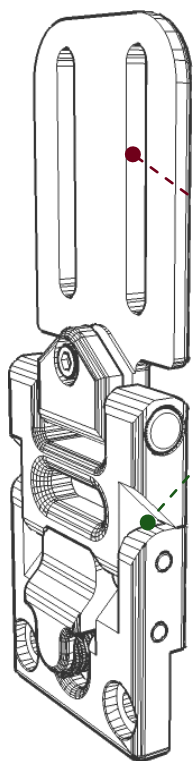

PAD MOUNT

Niño Lateral Pad Mount—TWBLNL-31

SWING-AWAY

Quick Release Swing Away Mechanism

QUICK RELEASE BUTTON

Spring Loaded Quick Release Button

TWBTL

Lateral Thoracic Thin Articulating Adjustable Swing Away Hardware (E1028)

The World's Best Thin Laterals are thin swing away thoracic lateral brackets that are designed with a thinner profile for mounting in areas where there isn't enough room for mounting the standard swing away lateral supports. The TWBTL comes with a non modular flip down bracket and a pad mount.

PRODUCT DESCRIPTION

Now the TWB Link hardware thickness has been reduced 25%! So the TWBTL Laterals 'snug' more easily and still retain all the strength, flexibility and function that make them the World's Best Laterals! With a variety of mounting options and modifications, Stealth's Thin Laterals are a versatile and strong while providing a low profile lateral for when space requires a thinner hardware to provide support.

TECHNICAL DATA

1. $\frac{3}{16}$ " tool for $\frac{8}{32}$ x $\frac{1}{4}$ socket head screws (hardware mounting.)
2. $\frac{5}{32}$ " tool for 10-32 x $\frac{3}{4}$ socket head screws (articulating joint.)


SELLING POINTS

1. Push button quick release swing-away mechanism provides easy transfer from user to chair.
2. Smooth finished edging on every piece.
3. Adjustable to accommodate the needs of almost any client.
4. Thin profile design for mounting in areas with little room for mounting adductors.

DIMENSIONS

Side View

Front View

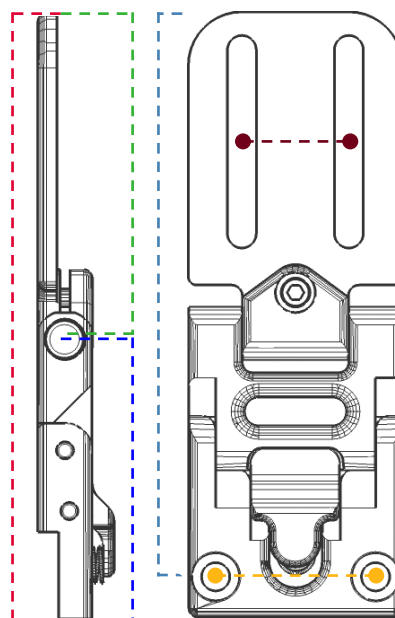


Standard Configuration Dimensions

Use These standard dimensions to determine which configuration you will need. See back for alternate mounting hardware and additional options.

Total Hardware Length	5½"
Pad Mount	3"
Swing Away Mechanism	2½"
Pad Mount Slot Spacing	1"
Swing Away Mechanism Hole Spacing	1½"
Mounting Hardware to Pad Mount	5¼"


2GTR Track Nut	LTHW503
Standard Tack Mount	LTHW515
All Season Bracket	LTHW504
Planar Mount	LTHW505
Formula Contour	LTHW507

HARDWARE FEATURES

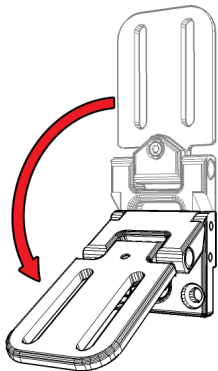

Swing-Away Action

Stealth's thin laterals include an easy to access button on the side of the hardware for the swing away mechanism. With a quick press of the side button the lateral swings out of the way to allow for easy user transfers.

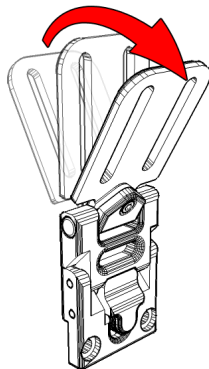

Adjustability

Thin links are extremely versatile hardware that are adjustable to provide excellent support and positioning. Screws connect the link's joints and allow for adjustment. Loosening the screws allows the links to be adjustable and tightening them makes the links become rigid to provide the best support possible. Stealth offers two different styles of thin laterals; one with a link and one without a link.


More information available at www.stealthproducts.com

MOUNTING

Mounting Hardware

Stealth offers mounting hardware of various sizes and style. Our diverse mounting hardware provides multiple options to mount on virtually any mobility base.

CANE MOUNT CLAMPS

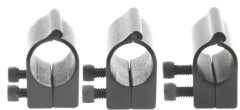

1" Cane Clamp—LTHW500

3/4" Cane Clamp—LTHW501

1/2" Cane Clamp—LTHW502

PLANAR MOUNTS*



All Season Bracket
LTHW504



Planar Mount
LTHW505

TRACKS



2" Track—SBHD-234

3" Track—SBHD-235

6" Track—SBHD-233

7" Track—SBHD-231

8" Track—SBHD-236

10" Track—SBHD-140

12" Track—SBHD-232

14" Track—SBHD-237

16" Track—SBHD-238

18" Track—SBHD-239



2nd Generation
Track Nut
LTHW503


RETROFITTING



Tru Comfort Back Mount
LTHW508



Formula Contour Bracket
LTHW507

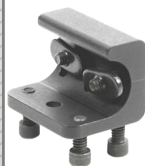


Pride Synergy Cane Clamp
LTHW506



Permobil Corups Tracks
Left—LTHW518-L
Right—LTHW518-R

ASSEMBLY DIAGRAM

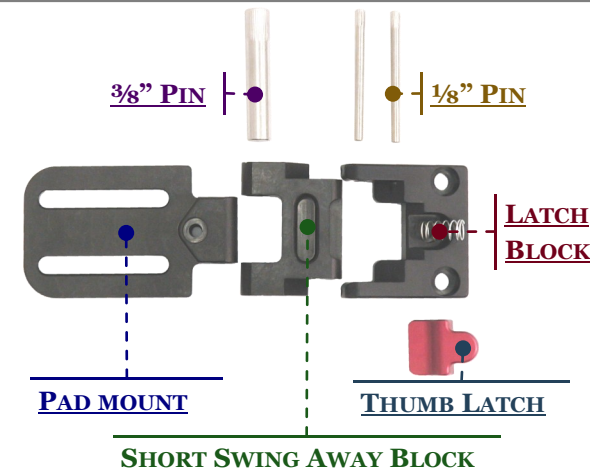

LATERAL PADS



TWBL3x5
3"x5" Niño Lateral Pad

TWBL4x6
2"x3" Niño Lateral Pad

Niño Lateral Pads

With Stealth Product's thoracic positioning/support pads the sky is the limit. We not only offer a wide selection of standard size pads made from an aluminum base with 1" closed cell foam. We can also make any other size one might dream up with foams and covers of choice.†

† When ordering different upholstery, foam or non-standard size the

*Both the standard and all season bracket are available for curved backs.